

NIRS DEVICES

X-RAY PHANTOMS

artinis

JOBS (/JOBS)

ABOUT US (/#ABOUT-US-SECTION)

NEWS (/#NEWS-SECTION)

PUBLICATIONS (/#PUBLICATIONS-HOMEPAGE-SECTION)


WE OFFER (/#WE-OFFER-SECTION)

CONTACT (/#CONTACT-HOMEPAGE-SECTION)

BLOG (/BLOG)

# Oxysoft

Our proprietary NIRS recording and analysis software.

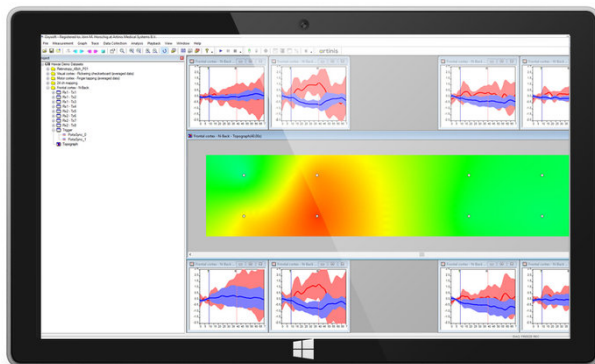

## Oxysoft

Provides in real time the calculation of oxy-, deoxy-, total hemoglobine and tissue saturation index (TSI)

- Applied to all of our NIRS devices
- Real time calculation of oxy-, deoxy-, total hemoglobine and tissue saturation index (TSI)
- Displays easy to understand 2D and 3D plots
- Supports simultaneous recordings of multiple devices
- Records and generates events from a wide range of sources
- All standard analysis techniques can be applied in realtime or posthoc
- Compatible with e.g. HOMER2, NIRS-SPM, FieldTrip, NIRXstorm, Polhemus and more
- Exports raw data in various formats
- Enables real-time streaming of data to e.g. Matlab or Python

Oxysoft is our proprietary NIRS recording and analysis software. Oxysoft allows robust recordings even of multiple devices at the same time, in any thinkable combination. This means you can combine recordings from several Oxymons, Portamons, Portalites, and Octamons within one measurement by using just one program. All data is synchronized and stored in one data file. You can analyze your data in real-time, during the recording, or also afterward.

## About Oxysoft

## Analyse your data

## All features


- It's main features include:
1. Data collection
  2. Import and export of Oxymon data files
  3. Exporting can be to a text file, Microsoft Excel or XML
  4. Import and scaling of data acquired via the optional analog input channels
  5. Importing of data files collected with other instruments
  6. Re-sampling of data
  7. Filtering of data using various types of filters
  8. Viewing of data
  9. Calculating average, mean, standard deviation, and regression
  10. Calculating differences between selected time-frames
  11. Calculation of oxygen consumption, blood flow, venous saturation, and arterial saturation
  12. Calculating averages of (block) stimuli (for functional NIR studies), including detrend function
  13. Generating events from external inputs, very useful for functional NIRS studies
  14. Preparing two-dimensional plots
  15. Preparing videos of the 2-D plots, e.g. for presentations. Very useful for brain mapping!
  16. Comparing between traces

20. Realtime data export to Matlab and other software

21. Offline data export to Matlab and Matlab based (SPM) software packages

22. Synchronisation to other software and devices

[SOFTWARE INFORMATION \(/S/OXYSOFT-LEAFLET.PDF\)](#)

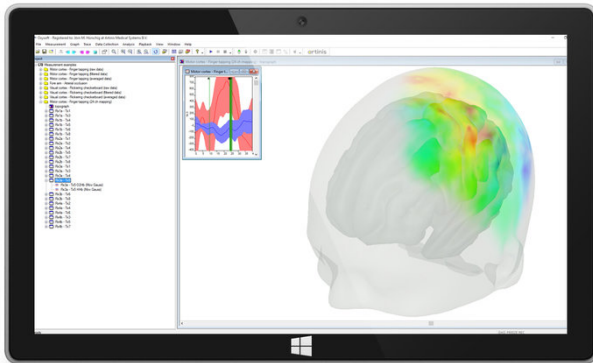

## OXYSOFT EXTENSIONS

We offer the opportunity to access extended visualization options, designed to optimize your measurement protocol, better review your results, and to optimally support your high impact neuroscience publication. The Oxysoft 3D extension can visualize your data on 3D models of the brain of the human head based on the widely used MNI templates.

With the 3D extension you are enabled to:

- Visualise your data on 3D models of the brain of the human head based on the widely used MNI templates.
- Use built-in, intuitive user-interface for the Polhemus fast-track device
- Automatically store the digitized positions of the optodes

[GET A QUOTE](#)

**OXYSOFT 3D & POLHEMUS**

CONTACT (/#CONTACT-HOMEPAGE-SECTION)

BLOG (/BLOG)

**Oxysoft 3D Coregistration and...**

We have partnered with Polhemus Inc, which creates devices for digitizing objects such as the Polhemus Patriot or Polhemus Fastrack device. Using these devices you can precisely measure the position of the optodes on the subject's head. Oxysoft 3D guides you through this process with a built-in, intuitive user-interface. No third-party software is needed for this process. Alternatively, you can import digitized positions from other manufacturers such as ANT's xensor or Magstim's EEG pinpoint device. Oxysoft automatically stores the digitized positions of the optodes, and transforms and then visualizes on the integrated 3D models. This not only allows you to better review your data and results but additionally allows to report back the optode positions in MNI coordinate space.

*Artinis NIRS products are intended to be used for research applications only. Our products are not sold as Medical Device as defined in EU directive 93/42/EEC. Our products are not designed or intended to be used for diagnosis or treatment of disease.*

## Be the first to know

Sign up for our newsletter and receive updates about the latest

NIRS NEWS

## Careers

We are hiring! Are you our new colleague?

news, publications and  
more.

NIRS DEVICES

X-RAY PHANTOMS

EVENTS

JOB

View all jobs

(/jobs)

ABOUT US (/#ABOUT-US-SECTION)

NEWS (/#NEWS-SECTION)

PUBLICATIONS (/#PUBLICATIONS-HOMEPAGE-SECTION)

WE OFFER (/#WE-OFFER-SECTION)

CONTACT (/#CONTACT-HOMEPAGE-SECTION)

BLOG (/BLOG)

CAREERS (/JOBS)

COMPANY

THE SCIENCE

BLOG (/BLOG)

SUPPORT (/SUPPORT)

CONTACT/ROUTE

Artinis Medical Systems, Einsteinweg 17, 6662 PW Elst, The  
Netherlands +31481350980 [hello@artinis.com](mailto:hello@artinis.com) (<mailto:hello@artinis.com>)

© 2000 - 2019 Artinis Medical Systems

([mail](mailto:hello@artinis.com)) ([http](http://www.artinis.com)) ([http](http://www.artinis.com)) ([http](http://www.artinis.com)) ([http](http://www.artinis.com))